
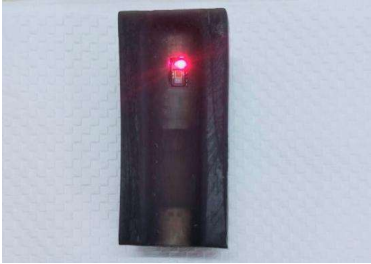
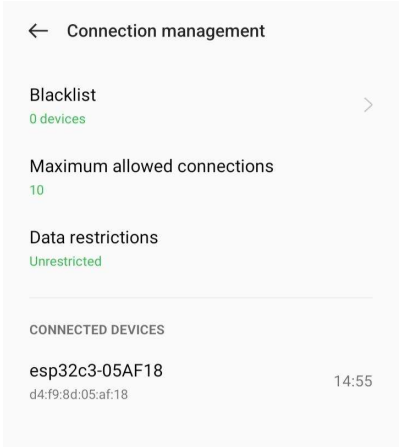


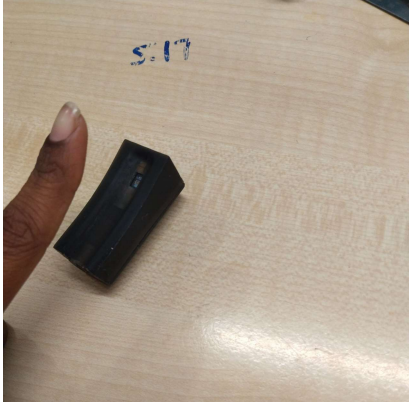
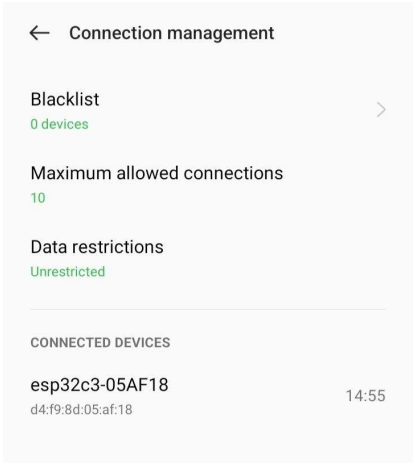


Sleep Apnea Monitoring Device Usage Protocol

S. NO	HEADING	IMAGE	CONTENT
1	Device Placement on the Finger	 <p>Device Off</p>  <p>Device On</p> 	<p>1.1 Preparation:</p> <ul style="list-style-type: none"> • Clean the patient's finger thoroughly to ensure good sensor contact. • Ensure the device is charged and powered on. • Red light is the indication for condition of the device. <p>1.2 Device connection confirmation</p> <ul style="list-style-type: none"> • Turn on the hotspot page while device is connecting to your mobile . • If the device is connecting to your mobile then , it will show the name of the device in your mobile. • After connecting the device , it will show like below this in your hotspot page for 1 to 5 seconds.

		<div data-bbox="448 210 866 544"></div> <div data-bbox="529 551 730 586"><p>Finger palcing</p></div> <div data-bbox="448 620 850 1014"></div> <div data-bbox="456 1021 847 1057"><p>Velcro Stricker Attachment</p></div>	<p>1.3Velcro Sticker Attachment:</p> <ul style="list-style-type: none">● Place the finger on the device properly .● Affix the Velcro sticker to the back of the device.● Wrap the Velcro strap around the patient's finger, securing the device in place.● Adjust the strap to achieve a snug but not tight fit to ensure comfort.
2	<p>Patient Comfort and Sleep Preparation</p>		<p>2.1 Comfort Guidelines:</p> <ul style="list-style-type: none">● Educate the patient on the importance of keeping the device securely fastened for accurate readings.● Emphasize that the device is designed to be non-intrusive to ensure a comfortable sleep experience. <p>2.2 Sleep Environment:</p> <ul style="list-style-type: none">● Advise the patient to create a sleep-friendly environment to minimize unnecessary movements during the test.

3	Handling Device Movement		<p>3.1 Repositioning Guidance:</p> <ul style="list-style-type: none"> • If the position is changed, the patient should check the device's position and reattach if necessary or contact the Technician for any help immediately. • Encourage the patient to be mindful of the device's placement, especially during sleep.
4	Finger Remove	 <p>Finger Removing</p>  <p>The screenshot shows a mobile application interface for 'Connection management'. It includes sections for 'Blacklist' (0 devices), 'Maximum allowed connections' (10), and 'Data restrictions' (Unrestricted). Below these is a 'CONNECTED DEVICES' section listing a device with ID 'esp32c3-05AF18' and MAC address 'd4:f9:8d:05:af:18', connected at 14:55.</p>	<ul style="list-style-type: none"> • If you want to remove your finger first take the velcro from the device. • While removing your finger then it will take the readings upto 1 minute. • If you are not placing the the finger for 1 minute , then it will wait for 3 minutes and after that the device is uploading our data into server. <ul style="list-style-type: none"> • While removing your finger from device, remember that it needs some mobile data to upload your data into ftp server . • Before removing your finger , you should turn on the hotspot in your mobile. • If the device is connected to your mobile , then it will be display the name of the device. • It is the indication for your data is uploaded into server. <p>Device Sleep:</p> <ul style="list-style-type: none"> • After uploading the data into server ,the device is going to sleep for 10 minutes. • After 10 minutes the loop is continuously executed.

5	Feedback and Alerts		<p>5. Feedback and Alerts:</p> <p>5.1 Device Status Indicators:</p> <ul style="list-style-type: none"> • The device may have LED indicators to convey its status. • Continuous On, Will Indicate the Device is in booting / Calibration. Blinking every 1 Second indicates proper connection and data collection, while continuous flashing may signal uploading the Data.
6	Emergency Situations:		<p>6.1 Device Disconnection:</p> <ul style="list-style-type: none"> • If the device disconnects due to significant movement or other issues, instructions for reconnection should be initiated by the Technician. • If reconnection is unsuccessful, the patient may need to contact technical support and / or reschedule the Test. <p>6.2 Unusual Patient Response:</p> <ul style="list-style-type: none"> • In the event of any unexpected discomfort, irritation, or adverse reactions, the patient should remove the device immediately and seek medical attention.

7.	. Patient Feedback and Reporting:		<p>7.1 Patient Diary:</p> <ul style="list-style-type: none"> Encourage the patient to keep a simple diary noting any significant events during the monitoring period, such as bathroom breaks or moments of wakefulness. <p>7.2 Patient Feedback Survey:</p> <ul style="list-style-type: none"> Provide a post-test survey to gather feedback on the patient's experience, comfort, and any challenges faced during the monitoring.
8	Follow-Up and Results:		<p>8.1 Post-Monitoring Consultation:</p> <ul style="list-style-type: none"> Schedule a follow-up consultation to discuss the monitoring experience, address any concerns, and explain the generated report. <p>8.2 Patient Education:</p> <ul style="list-style-type: none"> Reiterate the importance of following the monitoring protocol for accurate results and effective Sleep Apnea management.